

Innovation, safety and technology
bring you remarkable results.



11 and a half weeks after Cooltech® Define treatment.



11 and a half weeks after Cooltech® Define treatment.




14 and a half weeks after Cooltech® Define treatment.



11 and a half weeks after Cooltech® Define treatment.

Before & After
ONLY ONE SESSION

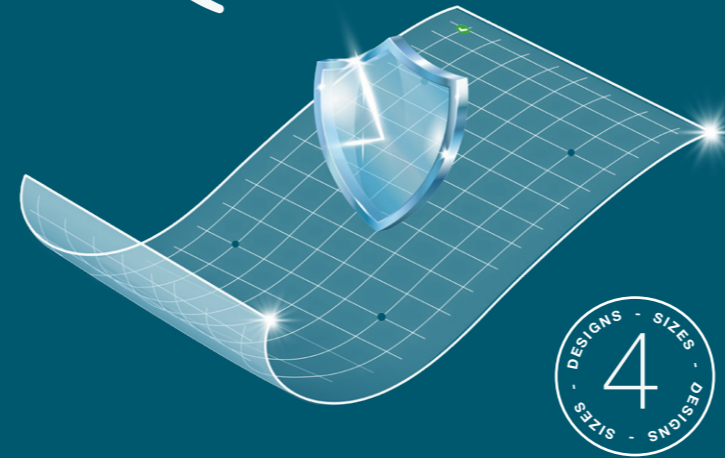

PATENTED
Cool Gel Pad
-50°C Smart Contact

The Cool Gel Pad -50°C Smart Contact is designed to protect the skin from the low temperatures during a cryolipolysis treatment.

New formula with -50°C gel which decreases freezing point and allows deeper penetration in the skin. Higher viscosity to keep the gel in the centre of the pad where the applicator will be placed protecting the skin during all the treatment.

Film located on the area in contact with the applicator (never in contact with the patient's skin) to avoid obstruction of the applicator by the gel and to facilitate applicator cleaning after treatment. Positioning markings for accurate placement of the applicator on the centre of the Cool Gel Pad -50°C Smart Contact.

TECHNICAL CHARACTERISTICS

Vacuum regulation level: 100 mbar - 500 mbar (± 10%)

POWER

Power: single-phase mains 100 - 240 VAC.

Mains frequency: 50 - 60 Hz.

Maximum power: 2400 W.

REGULATION

Classification according to IEC60601-1:


Class I, Type B (cooling applicators), Type BF (massager applicator)

PHYSICAL CHARACTERISTICS.

Dimensions: 715 x 1160 x 555 mm.

Weight: 90 kg.

Graphical interface: 10.4-inch touch screen

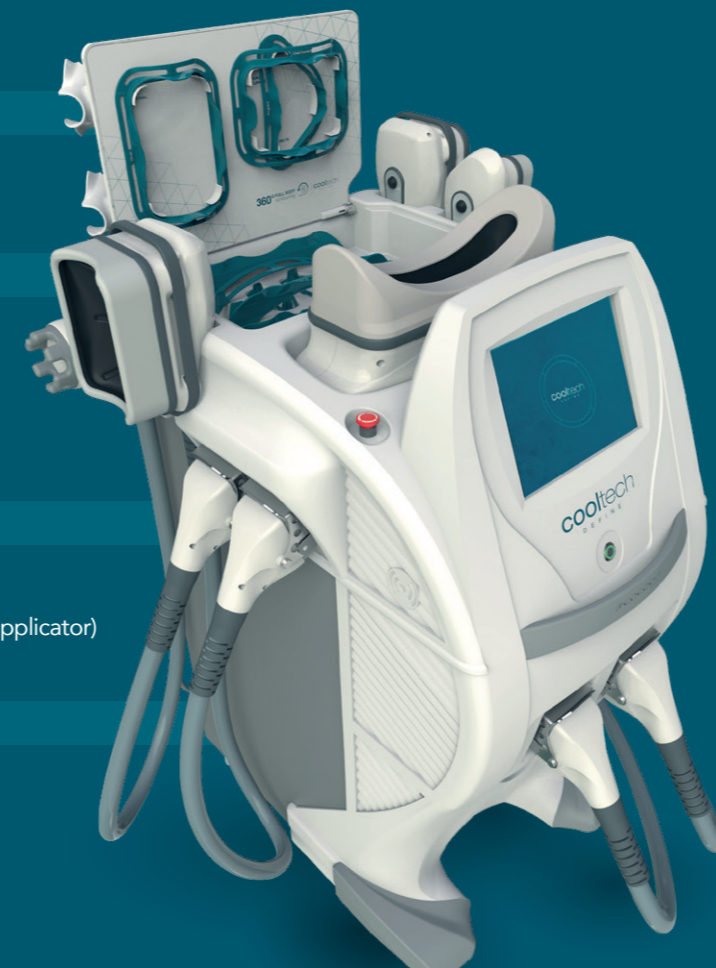

100014781 - 20230317 - V.5

cooltech®
DEFINE

FULL BODY
CONTOURING
LOCALISED FAT REDUCTION



Cooltech® Define is a cutting-edge cryolipolysis device that allows the treatment of **4 AREAS SIMULTANEOUSLY**

360°
COOLING
PLATES

Produces superb results, so every patient will want to experience this safe and highly effective treatment.

ebd@sinclair.com
+34 934 588 566
www.sinclair.com



OFFER COMPLETE BODY
CONTOURING TREATMENTS
IN RECORD TIME

**360° COOLING
PLATES**

Total contact of the plates with the tissue for better fat freezing results.

**9 ERGONOMIC
APPLICATORS**

Applicators designed to adapt to all areas, tissues and needs.

**4 APPLICATORS
WORKING IN
UNISON**

Treat 4 areas of your patients at the same time. Make the most of your time and their time.

**9 DIAGNOSTIC
TEMPLATES**

Ensure the best results with accurate assessments and optimise planning time.

**Patented Cool Gel Pad
-50°C Smart Contact**

Protect your patient's skin throughout the whole treatment.



MAKE THE MOST OF YOUR BUSINESS


LOCALISED FAT REDUCTION



TOTAL CONTACT
TOTAL COOLING
TOTAL RESULTS

**9 Ergonomic
applicators**

Tailor the treatment to all your patients.

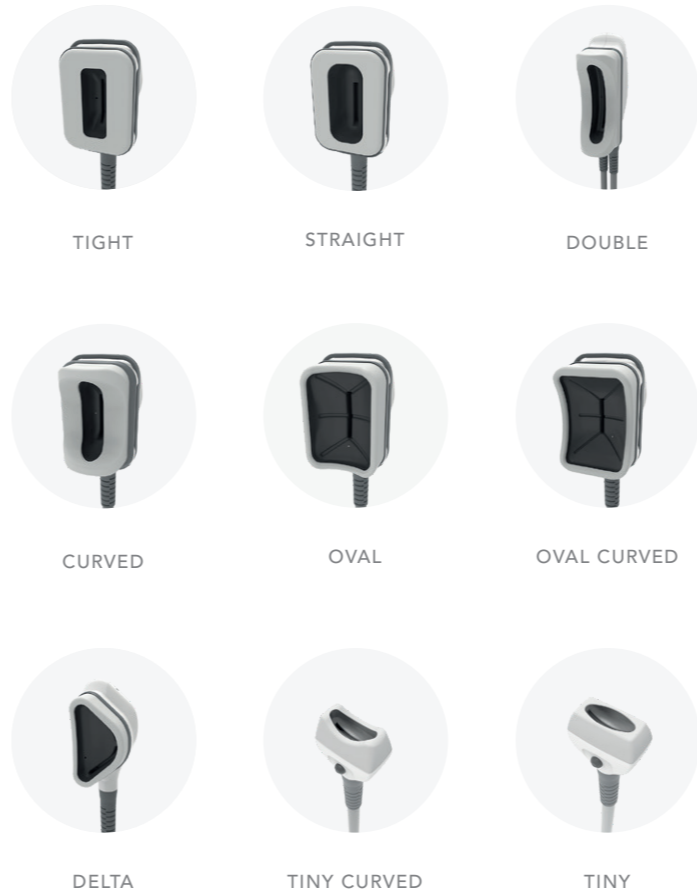


1 Massager

Boosts the efficacy of the treatment as it restores post-treatment blood circulation.




**Applicators
working
in unison**

Multiple combinations in just one session.

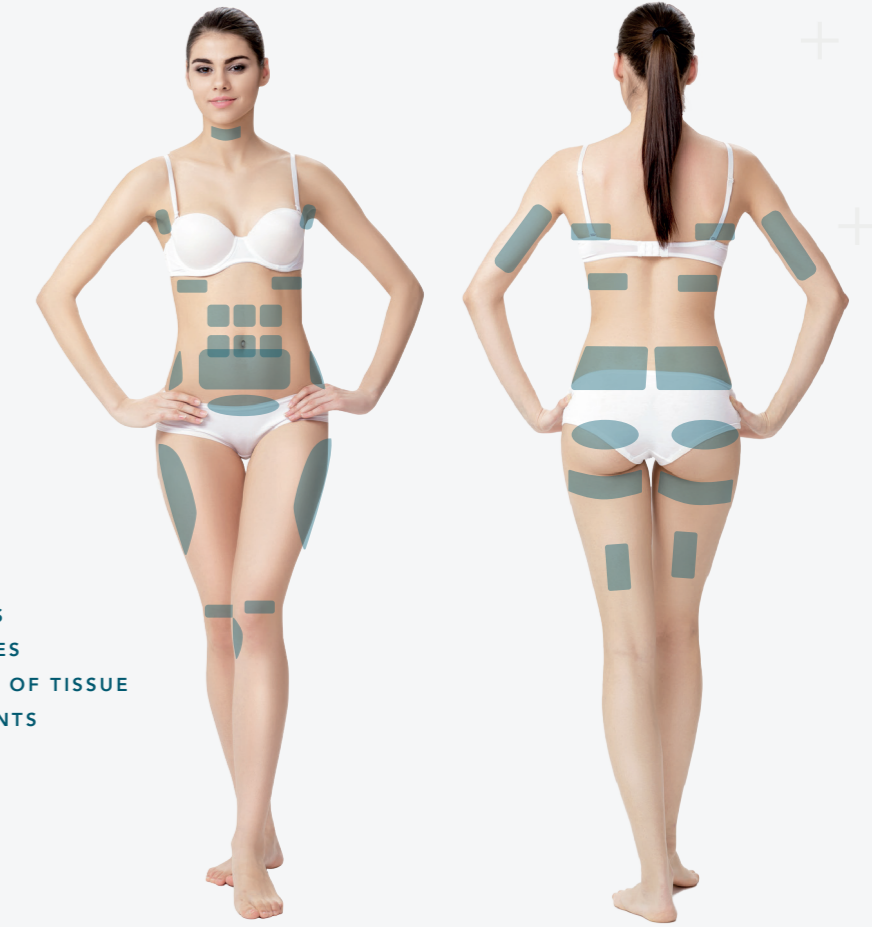

ALL AREAS
ALL ANGLES
ALL TYPES OF TISSUE
ALL PATIENTS

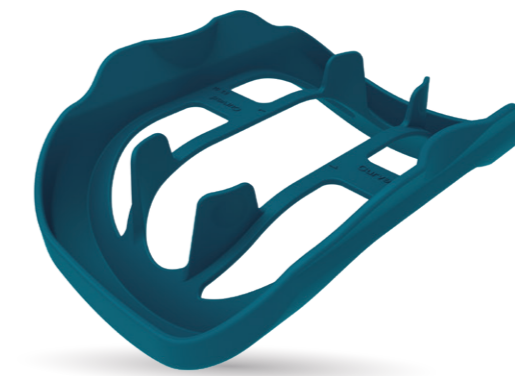

**User-friendly
diagnostic templates**

Ensure the best assessment and optimise planning time.

